

# Bobath XB Mat Table

Sturdy treatment table for neurological therapy



Main features:

1. Very high safe working load (SWL) 400 kg
2. Patented structure that makes the table sturdy and well suited for heavy use
3. Patient lifts fit easily underneath even when the table is at its lowest height (40 cm)

The renewed Lojer Bobath Mat Treatment Table is designed primarily for neurological therapy and Bobath therapy. The table has a patented structure that makes it very sturdy and well suited even for constant heavy use. It is designed in such a way that the therapist can be on the table at the same time as the patient who is being treated - safe working load (SWL) is up to astonishing 400 kg.


The tabletop is upholstered to a high standard and the sides are strengthened with plastic edging. The upholstery material is extremely durable and antimicrobial fabric that is also easy to clean, waterproof and has excellent disinfectant features. Tabletop's foam thickness is 30mm, which makes the patient surface less elastic and more suitable for examinations and treatments of neurological disorders. Bobath's frame structure allows patient lifts to fit easily underneath it (ground clearance is 15 cm) even when the table is at its lowest height. Extremely low minimum height (40 cm) ensures that even patients with wheelchair are able to get onto the table with ease.

There are two models available: 4010XB has one-part tabletop and 4020XB comes with two-part tabletop and adjustable back section. Both models can be equipped with optional face hole to enhance patient comfort.

**LOJER**® For easy care

## Technical data and features

| Model                                  | 4010XB   | 4020XB   |
|--|--|--|
| Sections                               | one-section  | two-section  |
| Total width                            | 1100 cm / 120 cm   | 100 cm / 120 cm  |
| Total length                           | 200 cm   | 200 cm   |
| Total length with face hole (optional) | 207 cm   | 207 cm   |
| Total weight                           | 145 kg   | 145 kg   |
| Height adjustment range                | 40 - 95 cm   | 40 - 95 cm   |
| Back section length                    | -  | 70 cm  |
| Back section adjustment range          | -  | 0-73°  |
| Trendelenburg / anti-Trendelenburg     | 20° / 20°  | 20° / 20°  |
| Castors                                | Ø 150 mm central locking twin  | Ø 150 mm central locking twin  |
| Other standard features                | Plastic edging to protect upholstery   hand control with safety lock functions (press two buttons for 3 sec) | Back section adjustment with hand control and safety lock functions (press two buttons for 3 sec)   plastic edging to protect upholstery |
| Upholstery                             | Several material and colour options, more information on request   | Several material and colour options, more information on request   |
| Safe working load (SWL)                | 400 kg   | 400 kg   |
| Ground clearance                       | 15 cm  | 15 cm  |
| Electric system                        | 24V  | 24V  |
| CE-marking                             | Yes  | Yes  |
| Warranty                               | 2 years  | 2 years  |
| Country of Manufacture                 | Finland  | Finland  |

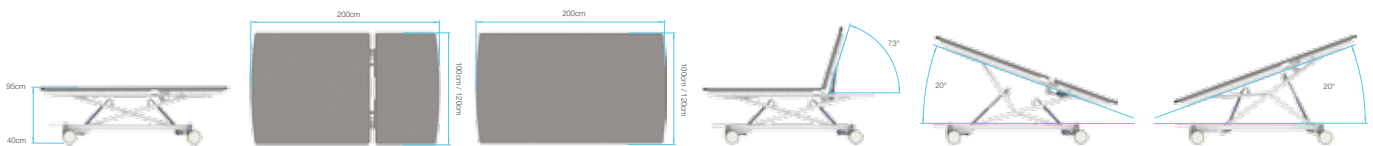

## Options

- Trendelenburg and anti-Trendelenburg (incl. battery)
- Foot control
- Face hole
- Battery powered

## Accessories

- Accessory rail
- Lift support
- Push handles
- Examination lamp
- Side rails



**LOJER**® For easy care