

LOJER®

*The human body features nature's high technology,
combining extremely complex functionalities
with unique personal properties.
This is why its repair and maintenance
is an expert task that requires special equipment.*

manuthera®₂₄₂

NEW DIMENSION IN MOBILISATION AND MANIPULATION TREATMENTS



More efficient,
accurate diagnosis

Excellent patient
comfort and safety

Improved efficiency,
high ergonomic standards

Features that enable
new treatment methods

Designed by experts, the most versatile treatment table in the world

Featuring unique modern technology designed by Finnish experts, the Manuthera® 242 mobilisation and manipulation table supports the therapist's work and ensures the best possible diagnosis and treatment for patients.

Manuthera 242 is an all-in-one treatment table for demanding examination and treatment use by various manual therapists, such as sports massage therapists, naprapaths, osteopaths, OMT and physical therapists and chiropractors. Manuthera 242 takes the therapist's work to a new level and supports countless treatment and diagnosis methods thanks to its patented operation system based on two synchronised motors, unique multi-traction features and various other functionalities, as well as its structure.


The patient can be moved to therapeutic and anatomical key positions smoothly and quietly. Mobility treatment can be carried out 3-dimensionally using the traction, flexion, lateral flexion and rotation features of the table, as well as gravity. The table design, twin motor system and the frame with bearings in all moving parts enhance the ergonomics of the therapist's work and make it less burdening, which enables even better treatment results.

Excellent treatment experience for the patient

Manuthera 242 offers an excellent treatment experience for the patient. Thanks to the sturdy frame, anatomically pleasant design and padding, as well as the easy-clean surface, the patient feels comfortable during all types of treatment.

The table movement is similar to that of the patient's spine, and there is no need for the therapist to lift or twist the patient. The smooth, easy movement of the table featuring natural routes enables painless, relaxing treatment. A pleasant treatment experience enables a higher number of repetitions and, consequently, better effect.

Thanks to the electronic tilting function, even patients with severe disabilities can easily get on the table. The technology is designed to support flexible adjustment to multi-method treatment and to make sure that the patient does not need to get off the table during treatment.

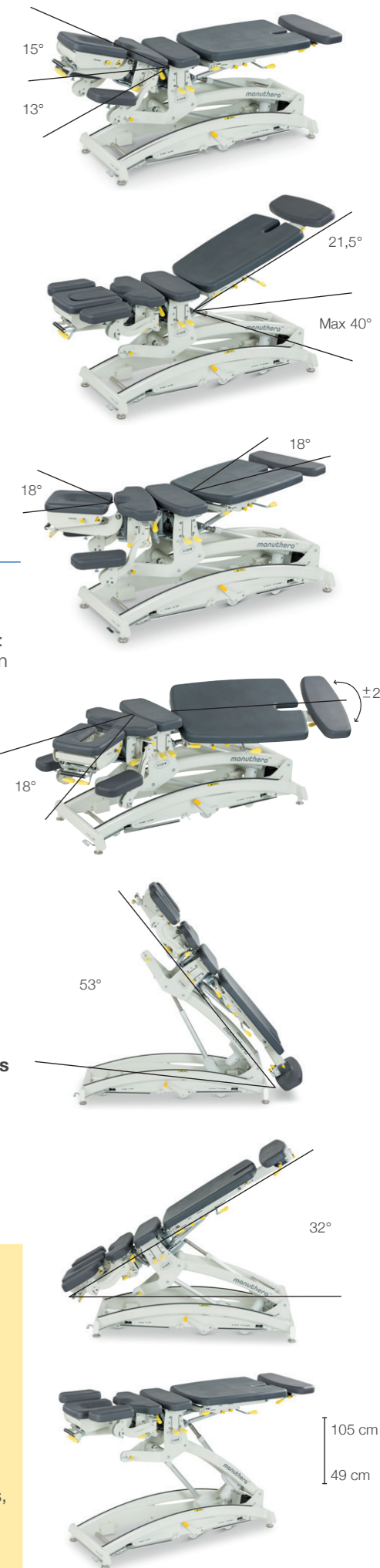

Supreme features for the best of the patient

Face section folding two ways and featuring advanced 3D functionalities for the safe mobilisation of the cervical spine and neck through vertical and lateral head movement.

Arm rests with stepless height adjustment and lateral movement. The arm rests can be attached to the side to function as side support.

Easy-to-use, sturdy and ergonomic **knobs** for locking the movement to the desired position.

Safety edge stops all table adjustments from being touched.

Adjustable **legs** with support paws.

NEW Two-piece thorax section element for vertical thorax bending and a 3D rotational motion. A unique tool for thorax mobilisation and manipulation!

Lower body section with extension/flexion, lateral movement and 3D rotational motion. **NEW** lower body rotation feature.

Multi-position foot section: upright position, footrest when lying down, and stepping board function.

Easy-care, durable, bio-compatible **soft moulded bed.**

NEW Drop function. New technology: traction drop, which enables easier and more diverse hip area manipulation.

Front and back **Total lock supports** provide additional support for burdening treatment techniques and make the table very sturdy in every section.

NEW Electronic Tilt up / Tilt down function for gravity-based treatment; also enabling easier access for patients onto the table. Spinal decompression has never been this efficient and pleasant for the patient!

Lockable castors for easy mobility.

Extremely wide range of **Hi-Lo adjustment.**

Manuthera 242 with synchronised lifting motors

Manuthera 242 is the world's first treatment table equipped with two synchronised lifting motors. This twin-motor feature enables smooth, easy movement and facilitates the therapist's work. The automatic safety feature detects if something becomes trapped between structures, which stops the motors and lifts the table to release the obstruction.

Manuthera® 242 – for the advanced professional

Manuthera 242 is designed to support the therapist's work in various ways:

- Timeless Finnish design for enhanced ergonomics
- Advanced treatment possibilities with multi-traction and 3D features
- Optimised motorised functions to reduce the need for burdening procedures
- Several movement options for each segment
- Functionalities to support treatment methods with a deep impact
- Ball bearings in table joints ensure smooth movement and endurance of high usage volumes
- Easy-care, durable surface materials
- Easy mobility of the table
- Pleasant treatment experience praised by patients
- High resale value
- Manufacturer's guarantee 2 years

	NEW	
	Manuthera 242	Manuthera 241
Thorax section elements	2	1
Height adjustment (Hi-Lo)	49–105 cm	53–89 cm
Weight	170 kg	128 kg
Tilt Down function (face section down)	Electronic	Mechanical
Tilt Up function (face section up)	Electronic	Mechanical
Lower body rotation	x	-
Total lock function for upper body	x	-
Total lock function for lower body	x	-
Castors	2	4
Safe working load (SWL)	150 kg	160 kg
Lifting motors (push only*)	2	1
Drop functions	Normal drop and traction drop	Traction drop
Automatic anti-squeeze safety function	x	-
Safety edge (movement stop)	x	-
Adjustable legs and paws	x	-
Safety lock functions	Automatic (2 h)	Manual with lever
Lying surface	Soft moulded	Foam plastic and artificial leather
U-shaped face element	x	-
Breathing hole cover	x	-
Country of manufacture	Finland	Finland


* Motor push to lift, gravitational lowering

Learn more about the product's features and value-added services, as well as the various financing options such as leasing, and find a retailer at: www.lojer.com


LOJER® For easy care

Россия, Москва, Вере́йская улица, 29с134
БЦ «Вере́йская Плаза 3» оф 306
Тел./факс:(495)646-80-28, почта: info@mtdiagnostica.ru

Manuthera 242
Maximum width 54 cm

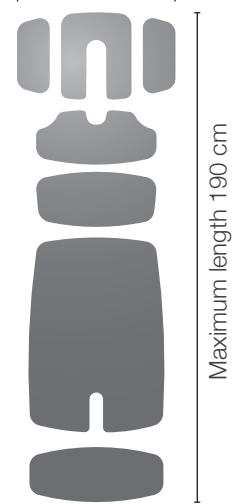

Manuthera 241
Maximum width 57 cm

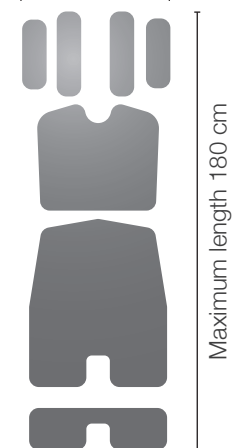

Also read about the Manuthera trainings in which the inventor of the product, naprapath Teuvo Keisala helps you make the best use of your table!



Manuthera is a registered trademark of Lojer Oy. Manuthera 242 mobilisation and manipulation tables are made by Lojer Oy, a Finnish company designing, manufacturing and marketing high-class tools and equipment for healthcare and wellness services in Finland and globally.